

snogg®

ANIMAL CARE



—Animal— WOUND CARE

FOR TREATMENT AND SUPPORT



Animal Wound Care since 1980

Snögg Animal Care has more than 40 years of experience in animal wound care. Developed in co-operation with veterinarians and animal wound care experts, the products are intended for animal use – with a focus on animal well-being. The product range comprises various types of bandages that are all easily adapted to different needs, areas of the body, and types of wounds, offering the user creative freedom.





Snögg Animal Care

- ✓ Easy to apply and remove
– without unnecessary stress for the animal
- ✓ Can be used to bandage difficult areas of the animal's body
- ✓ Multi-use – combine and customize for all animal types and stages of the wound healing process
- ✓ Can easily be combined with other wound care procedures
- ✓ For treatment and support
- ✓ Packaging design with clear instructions



Where can *Adhesive Foam Bandages* be used? Find out the answer and more on our website. Scan the QR code to visit our FAQ page.

Adhesive FOAM BANDAGE

Adhesive Foam Bandage is an elastic foam dressing with great adhesion to both fur and skin. The dressing has adhesion on one side and is particularly suitable for padding and fixation, even on areas that are difficult to bandage.

The product can be used on its own or in combination with a compress and/or **Thick Foam Bandage** to create a complete and stable wound dressing.



LATEX
FREE

SOFT AND FLEXIBLE

- Soft and flexible
- Adhesion on one side
- Suitable for padding and fixation

USAGE TIPS

- Fixation of compresses
- Fixation of ECG-pads
- Additional bandage fixation



*When you need great
adhesion to fur*

Adhesive Foam Bandage is available in three sizes to suit different animals and wounds:



REF 510002410
9 x 200 cm




REF 510002412
14 x 200 cm




REF 510002413
22 x 200 cm


BANDAGING WOUNDS:

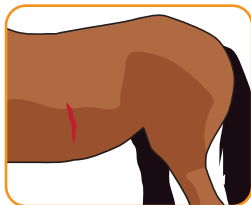

Ensure the wound is clean and the fur is dry.

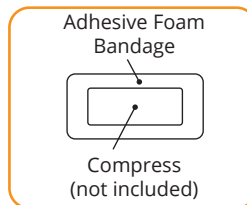


Use a suitable compress.

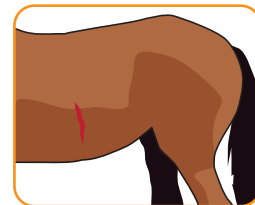

Fix the compress with the Adhesive Foam Bandage.


USING ADHESIVE FOAM BANDAGE TO PREVENT SLIPPING:

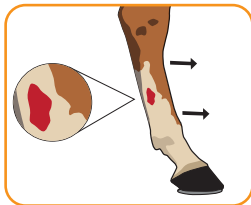

Apply Adhesive Foam Bandage below the wound before applying the compress and the outer bandage. The Adhesive Foam Bandage acts as a stop and helps prevent the bandage from slipping down, particularly on legs.



KNEE INJECTION:

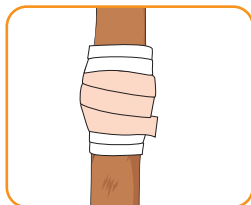

Begin by wrapping Adhesive Foam Bandage around the area below the knee. Then wrap Thick Foam Bandage over the injection point and overlap the Adhesive Foam Bandage. This acts as a stop and helps prevent the bandage from slipping down.

Additional fixation: For enhanced stability, Adhesive Foam Bandage can be applied over the bandage as an extra fixation layer. Especially suitable for highly mobile areas.

Scan the QR code for further instructions and instructional videos.



Thick FOAM BANDAGE

Thick Foam Bandage is a self-adherent foam dressing that only sticks to itself and therefore does not stick to wounds, fur or skin. The dressing is flexible, easy to shape and provides great mobility. The absorbing and soft foam makes Thick Foam Bandage particularly suitable for protection, padding and fixation.

REF 510002411
9 x 200 cm



**LATEX
FREE**



SOFT AND FLEXIBLE

- Self-adherent, only sticks to itself
- Can be placed directly on the wound
- Absorbing, soft and flexible


USAGE TIPS

- Protection of superficial wounds
- Paw bandaging
- Bandage highly mobile areas
- Protection during mould treatment




*Flexibility and
freedom to move
with Thick Foam
Bandage*


PAW BANDAGING:

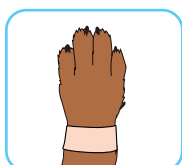


Wrap one length of Thick or Thin Foam Bandage around the leg.

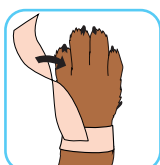

Tear off a suitable length to fit between the digits. Fix it to the wrap, turn it 180 degrees.

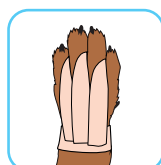


Slip the narrow part between the digits and fix the wrap on the other side. Repeat for each digit.

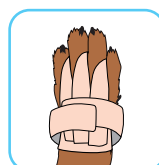


Complete by wrapping a second bandage around the leg to secure the bandage.


LEG BANDAGING:

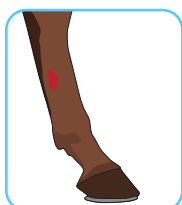

Ensure the wound is clean and the fur is dry.

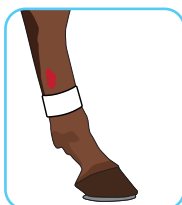

Begin by wrapping one layer of Adhesive Foam Bandage below the wound area for fixation.

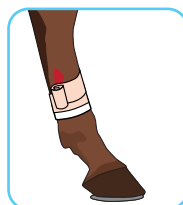

Then wrap Thick Foam Bandage around the wound in a spiral pattern. Overlap the first bandage by 2-3 cm.



Replace the second bandage when required. The first bandage is removed when healing is completed.

EAR BANDAGING:



Shape a length of Thick Foam Bandage into a cone and place it in the ear to protect the ear canal.



Fold another piece to create a triangle and place behind the ear on the intended wound or to support the ear when bandaging.




Cut the Thick Foam Bandage to a 4.5 cm strip and wrap around.



Bandage the ear with Thick Foam Bandage. Hold it in place with a wider Adhesive Foam Bandage that sticks to the skin and fur. Leave an opening for the ear canal and leave the healthy ear free of bandage.

MOULD TREATMENT:

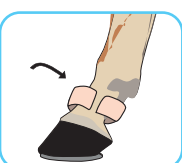

Thick Foam Bandage can be used as a protective layer during mould treatment. Mould requires regular monitoring.

Scan the QR code for further instructions and instructional videos.



Thin FOAM BANDAGE

LATEX
FREE



Thin Foam Bandage is flexible and easy to shape. The cross-elastic bandage sticks to itself and does not stick to the skin, fur or wound. It is soft, absorbent and can be placed directly on the wound.

REF 510002379
6 x 450 cm


SOFT AND FLEXIBLE

- Self-adherent, only sticks to itself
- Can be placed directly on the wound
- Absorbing, soft and flexible


USAGE TIPS

- Paw bandaging
- Teat bandaging
- Tail bandaging
- Protection of minor superficial wounds


PAW BANDAGING:

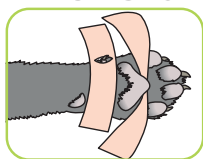


Place a length of bandage between the digits and the pad. Make a slot in the bandage and place it around the dew claw.

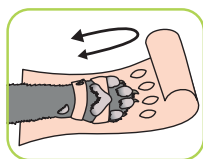

Tear off a suitable length of bandage. Make slots to go over the digits.

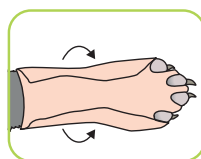

Slip the digits through the slots and fold the bandage as shown.


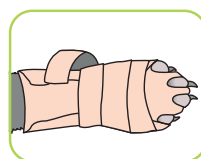


If necessary, wrap an additional length of bandage around the paw to secure the bandage.



Scan the QR code for further instructions and instructional videos.

TEAT BANDAGING:

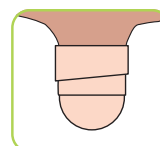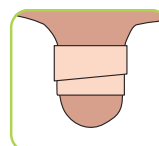

Use Thin Foam Bandage in the event of injury of udders or fixation of teat plugs.





*Elastic and self-adherent
Thin Foam Bandage*

Primary layers

Snøgg Animal Care offers two types of dressings for different kinds of wounds, including small wounds and heavily exuding wounds. They work very well when combined with other wound care products and are suitable for fixing with **Adhesive Foam Bandage**.


Foam COMPRESS



PROMOTES WOUND HEALING

A sterile and absorbent foam compress that keeps the wound moist and promotes wound healing. A soft and flexible dressing that follows the body's movement. It does not stick to the wound and can be cut to the desired size.

- Sterile and absorbent
- Provides moist wound healing
- Soft and flexible

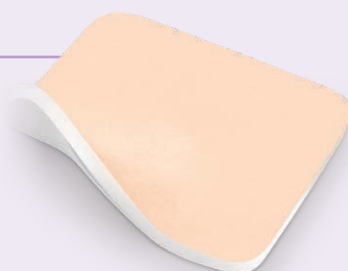

REF 510003122
10 pieces (10 x 10 cm)

Superabsorbent COMPRESS



FOR HEAVILY EXUDING WOUNDS

A sterile compress with high absorption capacity for heavily exuding wounds. The internal superabsorbent core can absorb a large amount of exudate. The exudate remains locked within the dressing pad to help prevent exudate from leaking through the dressing.

- Sterile
- High absorption capacity
- Helps to prevent leaking of exudate



REF 510003124
10 pieces (10 x 10 cm)

Product overview



	ADHESIVE FOAM BANDAGE	THICK FOAM BANDAGE	THIN FOAM BANDAGE
Sticks to itself – not to fur, skin or wounds	✗	✓	✓
Excellent adhesion to fur & skin	✓	✗	✗
Cross-elastic	✓	✓	✓
Length-elastic	✗	✓	✓
Absorbent	✗	✓	✓
Stops bleeding	✗	✓	✓
Adherent even when wet	✗	✓	✓
Suitable for primary layers	✗	✓	✓
Suitable for secondary layers	✓	✓	✓
Suitable for tertiary layers	✓	✗	✗
Fixation	✓	✓	✓
Thickness	5 mm	5 mm	1,9 mm



	FOAM COMPRESS	SUPERABSORBENT COMPRESS
Low exuding wounds	✓	✗
Medium exuding wounds	✓	✓
Heavily exuding wounds	✗	✓
Surgical wounds	✓	✓
Moist wound healing	✓	✗
Sterile	✓	✓
Can be cut to desired size	✓	✗
Suitable for fixation with Adhesive Foam Bandages	✓	✓



ANIMAL CARE

Distributor:

www.snogg.com

Orkla Wound Care AB (Svetsarvägen 15), Box 1336, SE-171 26 Solna, Sweden
Tel +46 (0)8 590 96325, snogganimal@orkla.com, www.snogg.com